



Creating and Editing Qualifier

To add a qualifier to an annotation select it in one of the *Sequence View* subviews and press the Insert key, or use the *Add Qualifier* context menu or the *Actions* main menu item.

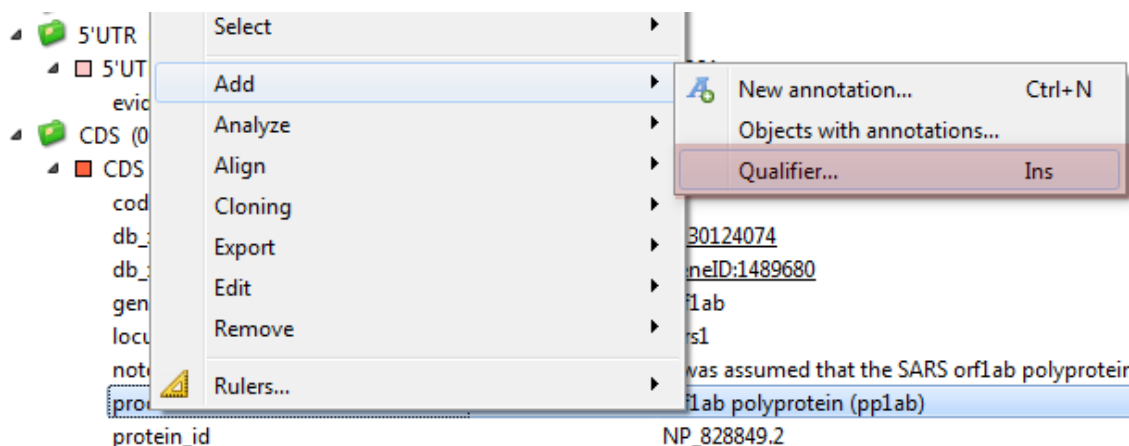


The dialog will appear:

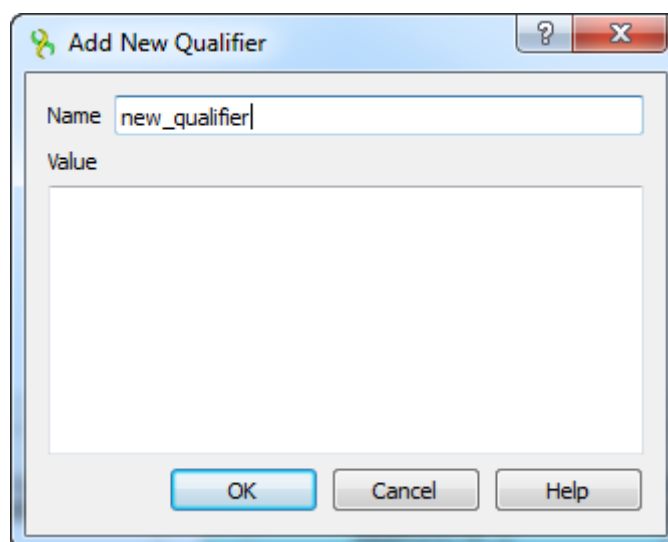

Here you can specify the name and the value of the qualifier.

You can use the F2 key to rename a qualifier:

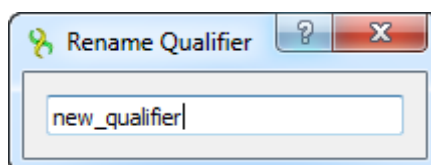

To edit a qualifier, select the qualifier and press the F4 key or use the *Edit qualifier* context menu item:

The image shows a software interface with a tree view on the left and a context menu on the right. The tree view is titled "Name" and "Value". It contains the following items:

- NC_001363 features [murine_copy2.gb]
- CDS (0, 4)
- misc_feature (0, 2)
- misc_feature (Value: 2..59)
- mi (Value: 45.)
- source

The context menu is open over the "misc_feature" item and contains the following options:

- Copy qualifier 'note' value
- Add 'note' column
- Edit qualifier** (highlighted with a red circle, with keyboard shortcut F4)
- Go to position... (Ctrl+G)
- Select sequence region... (Ctrl+A)
- New annotation... (Ctrl+N)
- Copy
- Select
- Add
- Analyze
- Align

At the bottom left, there is a "Tasks" panel with a magnifying glass icon and the text "2: Tasks".