



# Export Sequences Associated with Annotation

In UGENE you can export a sequence associated with an annotation. To do it select the annotation in the *Project View* window and click the *Export/Import* *Export corresponding sequence* context menu item:

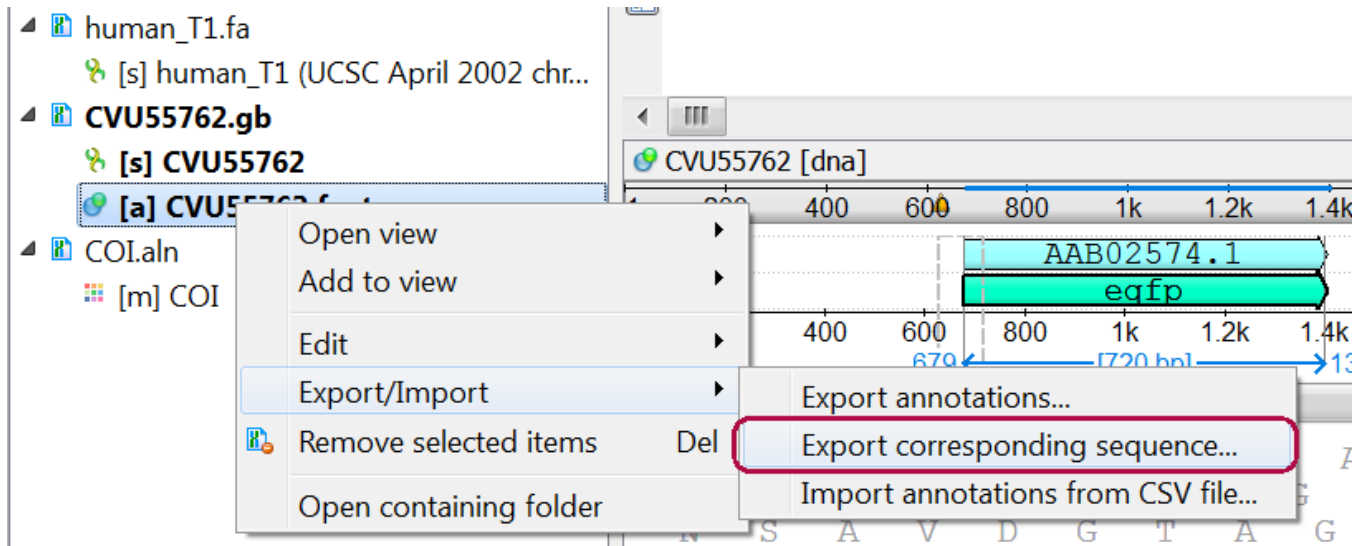

The *Export Selected Sequences* dialog will appear:

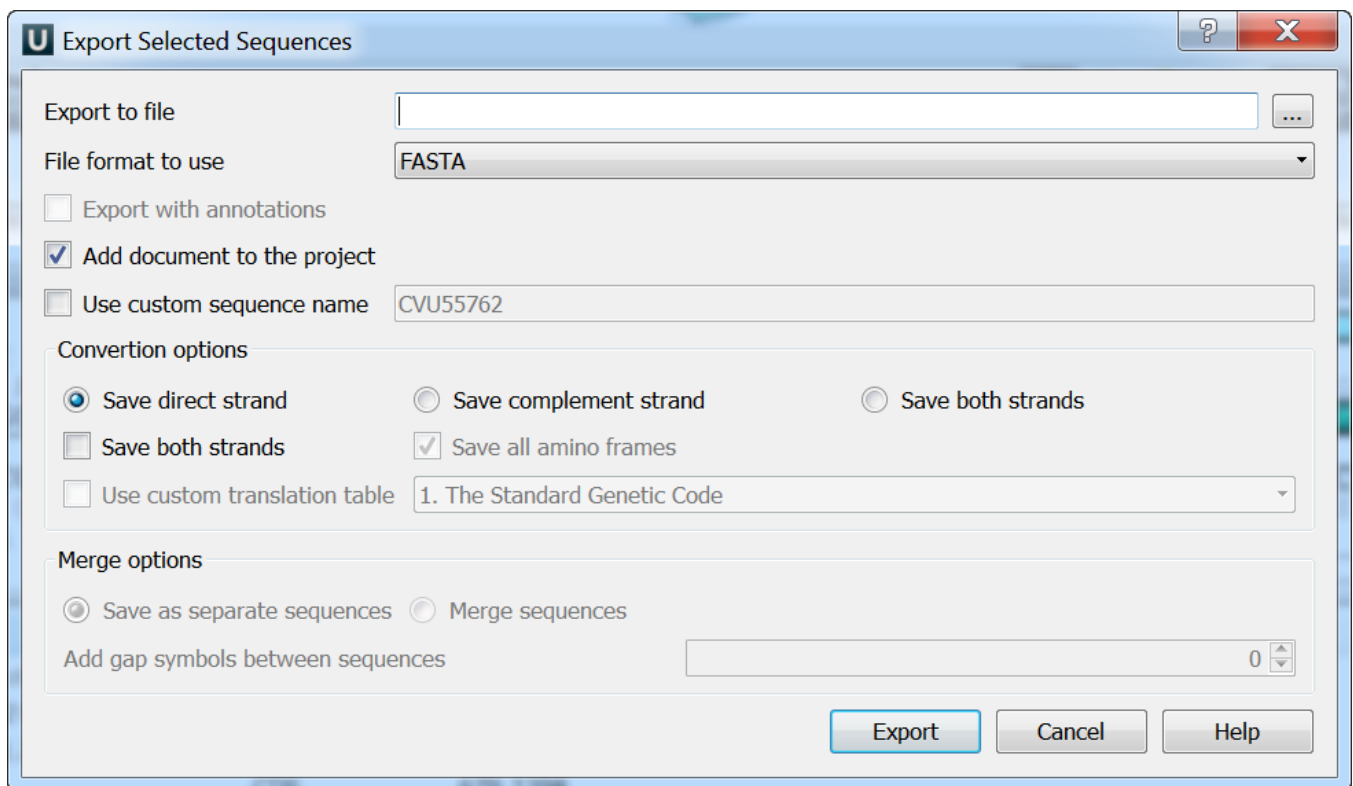

Here you can select the location of the result file and a sequence file format. You can choose to add newly created document to the current project and use custom sequence name. To do it check the corresponding checkboxes.

Use the *Conversion options* to choose a strand for saving sequence(s). Also you can translate sequence(s) to amino alphabet.

Also it is possible to specify whether to merge the exported sequences into a single sequence or store them as separate sequences. If you merge the sequences, you're allowed to select the gap symbols between sequences. This is the length of the insertion region between sequences that contain **N** symbols for nucleic or **X** for protein sequences.

## Export sequence with annotations

To export sequence with annotations choose Genbank or GFF format. The *Export with annotations* checkbox will be available. Check the checkbox and sequence will be exported with annotations .