


Multiple Sequence Opening

To open several sequences use the *File->Open* menu item or *Open* toolbar button and using *Ctrl* select the several sequences and click the *Open* button. The following dialog will appear:

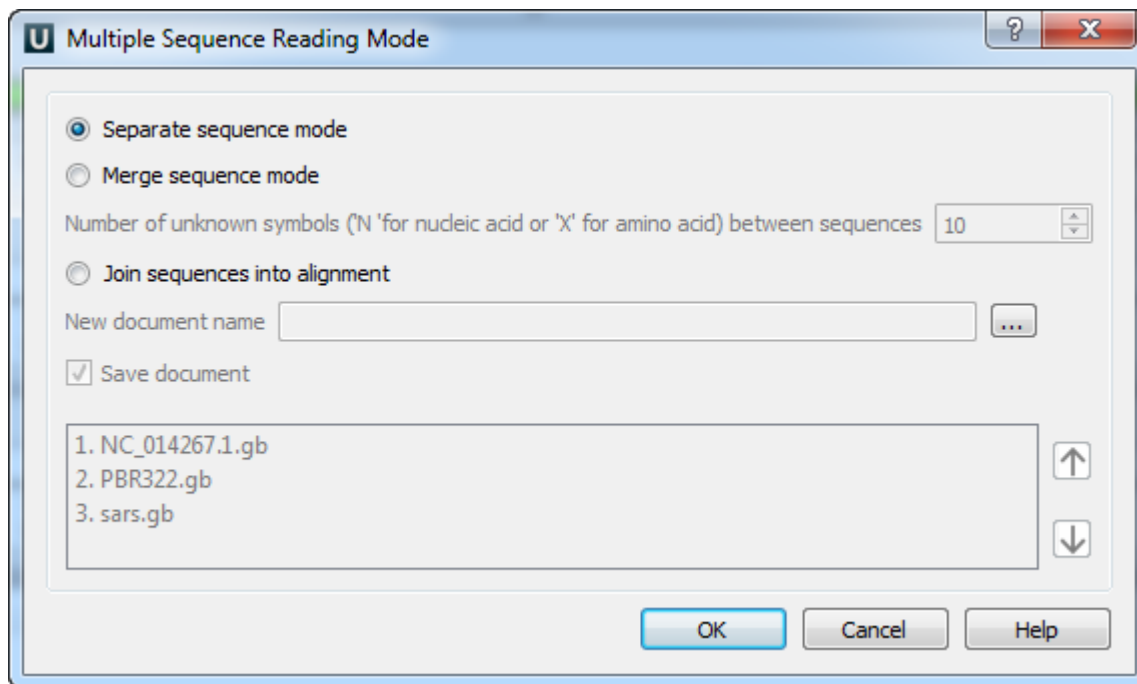

The following parameters are available:

Separate sequence mode - opens the sequences as separate sequences.

Merge sequence mode - merges sequences into one sequence with selected number of unknown symbols between sequences.

Join sequences into alignment - joins sequences into alignment.

Save document - save document to the selected document.

Also you can change the order of the sequences by up and down arrows.

Choose the parameters and click the *Open* button. Sequences will be opened in the selected mode.

In the *Separate sequence mode* sequences will be opened as separate sequences in selected order. You can change the sequences order by drag and drop in the sequence view.