


Export Chromatogram Consensus

To export consensus sequence use the *Consensus* tab of the *Options Panel*:

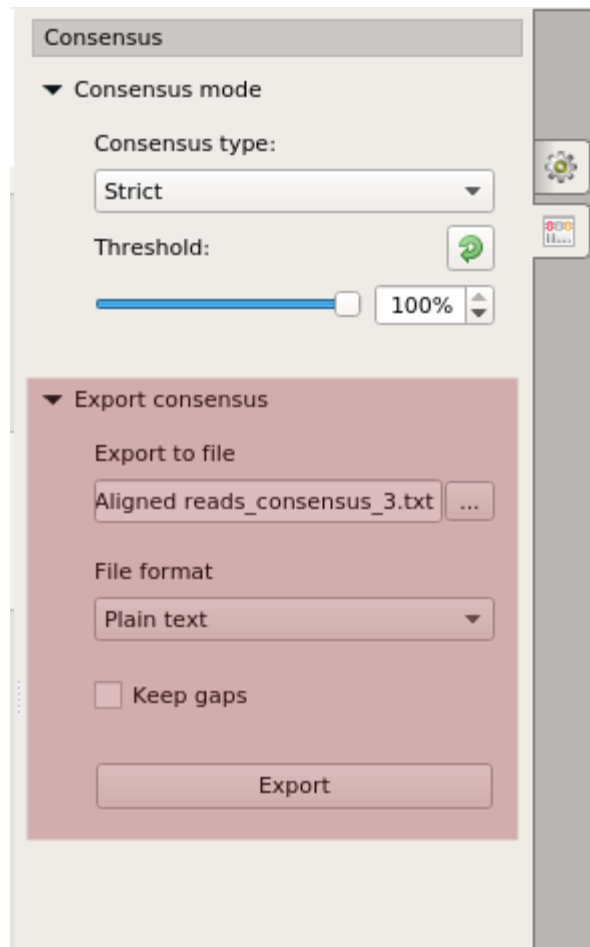

The image shows a software interface for exporting a consensus sequence. It is titled "Consensus" and has a dropdown menu for "Consensus mode". Under "Consensus mode", there are two settings: "Consensus type" is set to "Strict" in a dropdown menu, and "Threshold" is set to "100%" with a slider and a refresh icon. Below this is a section titled "Export consensus" with a red background. It contains "Export to file" with a text input field containing "Aligned reads_consensus_3.txt" and a file selection icon "...". Below that is "File format" set to "Plain text" in a dropdown menu. There is also a checkbox labeled "Keep gaps" which is currently unchecked. At the bottom of this section is an "Export" button. On the right side of the panel, there are two icons: a gear icon and a "BIO" logo.

The following parameters are available:

Export to file - here you need to select path for the output file.

File format - format for the output file.

Keep gaps - set this parameter if the result consensus must keep gaps.

When you click on the *Export* button the consensus sequence will be exported into selected output file.